

[home](#) [client](#) [applet](#) [daemon](#) [about](#)

## Users' Manual for Talairach Client 2.4.3

### Table of Contents

1. [About Talairach Client](#)
2. [Search Options](#)
3. [Searching](#)
4. [Command Lines Options](#)
5. [Suggested Uses](#)

Last Updated March 2009 (CLI gray matter search)

## About Talairach Client

The Talairach Client is used to assign Talairach Atlas labels for a given x,y,z coordinate. The assigned label is hierarchical, and is composed of five levels: hemisphere, lobe, gyrus, tissue type, and cell type.

### For example:

Right Cerebrum, Temporal Lobe, Sub-Gyrus, Gray Matter, Hippocampus

Inter-Hemispheric, \*, \*, White Matter, Corpus Callosum

Left Cerebrum, Temporal Lobe, Inferior Temporal Gyrus, Gray Matter, Brodmann area 20

A list of Talairach labels segregated by each level of the hierarchy can be seen [here](#), or a list of the entire set of labels can be seen [here](#).

## Search Options

A label search offers three options: single point, nearest gray matter, and cube range.

Database Search Options

☐ Single Point
☐ Nearest Gray Matter
☒ Cube Range

+/- 2 mm

### Single Point Search

A single point search is the simplest; it returns the label assigned to the given x,y,z coordinate.

### For example:

Label at (-20, 15, 5):

Left Cerebrum, Sub-lobar, Lentiform Nucleus, Gray Matter, Putamen

### Cube Range Search

A cube range search returns all labels within a cube centered on the given coordinate. The width of the cube is defined by the user, ranging from "+/- 1mm" (3mm wide) to "+/- 5mm" (11mm wide). Multiple voxels within the cube are likely to have the same label. The number of "hits" that each label gets is reported before each label. The sum of all of the hits is the volume of the chosen cube.

### For example:

Labels within 2 mm of (45, 45, 15):

84: Right Cerebrum, Frontal Lobe, Middle Frontal Gyrus, White Matter, \*

[Download Mac](#)

[Download PC](#)

[Download Other](#)

### Talairach Project Funding

The Talairach project has received support from the [EJLB Foundation](#) and the [NIH Human Brain Project](#).

This work forms part of a continuing project of the [ICBM](#) to develop a probabilistic atlas of human neuroanatomy.

### Recommended System Requirements

- [Java 1.4](#) or better
- 256 MB RAM

### Java Runtime Environment

Update the [Java Runtime Environment \(JRE\)](#) at Sun's official site.

### Feedback

Please email [Mick Fox](#) if you have any questions or find any bugs while using this software.

20: Right Cerebrum, Frontal Lobe, Middle Frontal Gyrus, Gray Matter, Brodmann area 46

16: Right Cerebrum, Frontal Lobe, Middle Frontal Gyrus, Gray Matter, Brodmann area 10


4: Right Cerebrum, Frontal Lobe, Middle Frontal Gyrus, \*, \*

1: Right Cerebrum, Frontal Lobe, Inferior Frontal Gyrus, Gray Matter, Brodmann area 46

### Nearest Gray Matter Search


Nearest gray searches involve concentric cube searches with varying diameters. In general, it searches consecutively larger cubes until it finds a gray matter label. It has the same outer limit of a 11mm wide cube, so it is possible to find no gray matter labels. **If two gray matter labels are found in the same search level, the Talairach Client returns the one with more hits. In the case of a tie, the next cube out is checked until the tie is broken.**

Since the gray matter search is limited to 11mm wide, in some cases it won't find a label. The result is then "No Gray Matter Found". If you would prefer to have a non-gray matter label instead of the "none found" message, choose that option in your preferences:

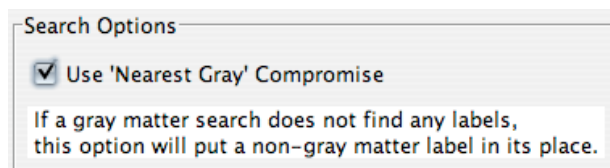


### Searching

To search for Talairach labels, you'll need coordinates. You can define the coordinates manually in the interface:

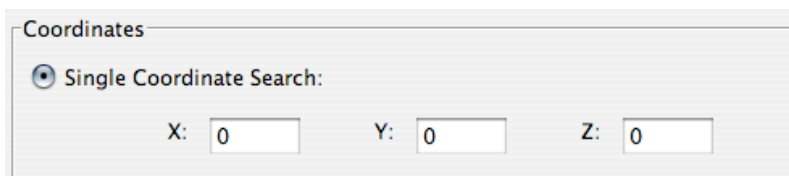


Or you can load coordinate data in from a file:

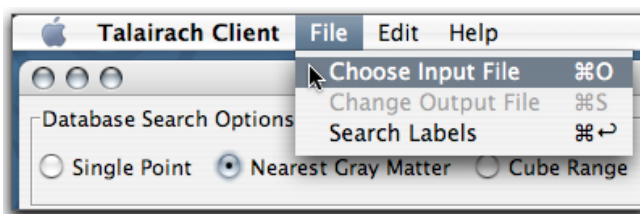


### Input Files

The Talairach Client can read text files delimited by tabs, spaces or commas. Multiple files can be chosen at the same time. Holding down shift will select multiple continuous files. To edit the selection more precisely, click while holding Control (on Windows) or Command (on Mac). Selecting multiple input files will create multiple output files, each one in the same directory with a generated name. For example, example\_data.txt became example\_data.td.txt. We put the ".td" before the extension so text and spreadsheet programs will recognize the file.

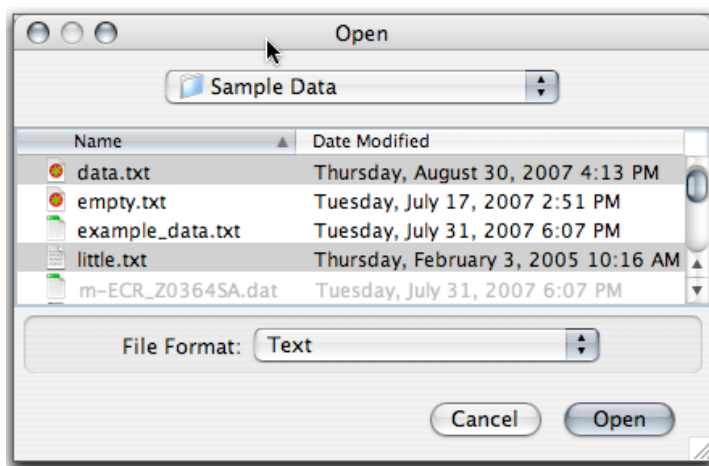


The Talairach Client also allows for header rows or extra columns before the coordinate data. To skip over extra rows or columns, you need to tell the Talairach Client which row & column is the start of your coordinate data. To assist with that, there is now an "Input Preview" window:

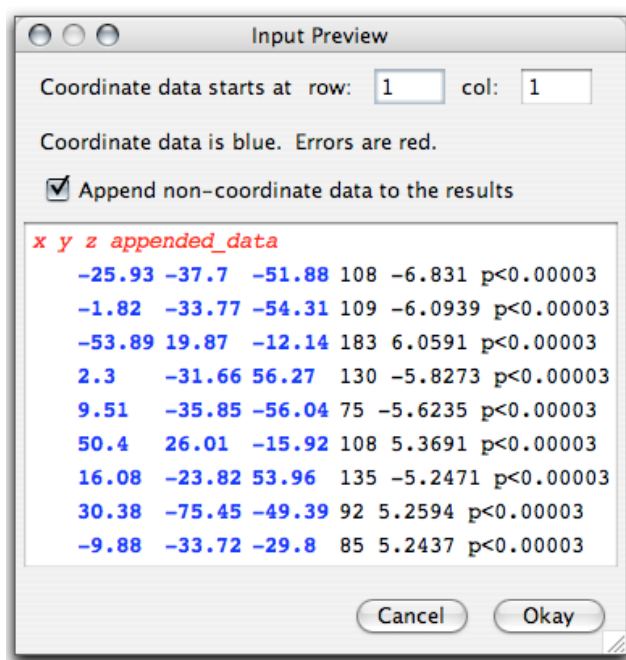

The Input Preview shows you the first ten lines of the file you have chosen. It tries to parse out the coordinate data, given the row and column number. If it has any problems, it displays the line in red. Coordinate data is highlighted with blue. There is also an option to "Append non-coordinate data to the results". If unchecked, only the coordinate data will be extracted. Otherwise, all data before and after will be appended to the results.

If you regularly open the same types of files with the same starting row and column, you probably will want to define a file type. The file type definitions are listed in the Preferences. The options are the same as those in the Input Preview with the addition of "Description" and "Extension".

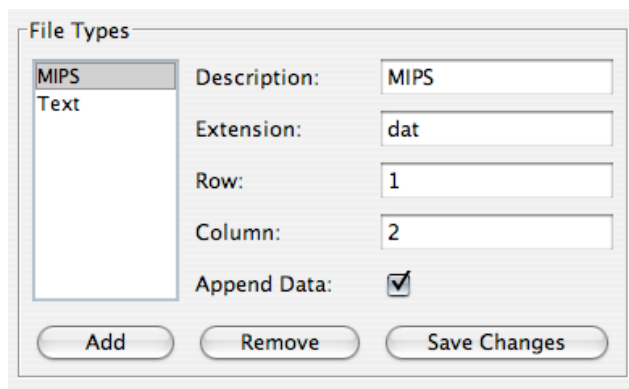


The "File Types" dialog box contains a list on the left with "MIPS" selected and "Text" below it. To the right, there are four input fields: "Description:" with "MIPS", "Extension:" with "dat", "Row:" with "1", and "Column:" with "2". Below these is a checkbox for "Append Data:" which is checked. At the bottom are three buttons: "Add", "Remove", and "Save Changes".

The description field is the name of this file type and how it will be listed in the "File Format" drop-down list when opening a file. The extension field lets you set a file extension. Only files that end with the chosen extension will be shown. Leave the extension field blank to allow all files. To define multiple allowed extensions, list them separated by spaces.

### Output Files

The standard format for Talairach Client output files is: record number, coordinate data, labels, and appended data.


The record number indexes each coordinate. The coordinate data from the original file is listed in the next three columns. Next, the Talairach label data is split up over several columns and changes depending on the chosen type of search. For a single label search, it is distributed over five columns. For nearest gray matter, it's the same five columns followed by a "range" column. Cube range reports "hits" first, then the five label columns. After the label data is any appended data from your original input file.

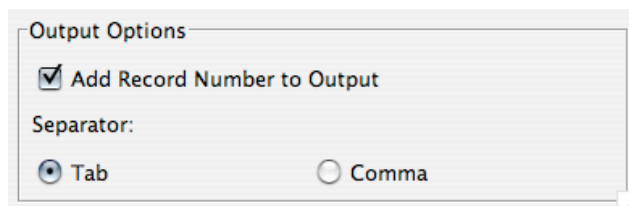

The "Output Options" dialog box has a checked checkbox for "Add Record Number to Output". Below it is a "Separator:" label with two radio buttons: "Tab" (selected) and "Comma".

Record numbers are optional and can be turned off in the "Output Options" section of the Preferences. Also, you can choose if the output files are tab- or comma-delimited.

If you are working with input files, the results will be saved to an output file in the same directory. You can view the results immediately you can click the "Show Results" checkbox. If checked, all of the text in the output file will be shown in the Talairach Client's text area.

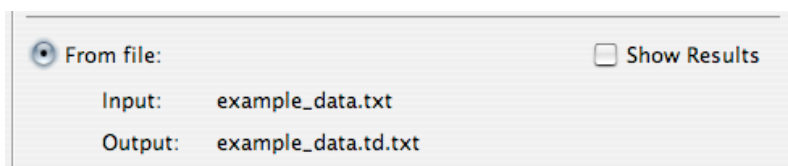

This section shows a radio button selected for "From file:" and an unchecked checkbox for "Show Results". Below, the "Input:" field is set to "example\_data.txt" and the "Output:" field is set to "example\_data.td.txt".

### Command Line Options

The Talairach Client includes two command-line tools, PointToTD and ExcelToTD, that communicate directly with the Talairach Daemon. By default, an internet connection is **required** to use these tools because they send queries to the Daemon. The host is talairach.org and the port used is 1600. If you have a [local Daemon](#), an internet connection is not required.

To use the command-line tools, you will need the [.jar](#) file.

## PointToTD

PointToTD operates on a single coordinate and has three search options: (1) Structural Probability Maps, (2) Talairach label, (3) Talairach labels within a cube range and (4) Talairach gray matter labels. The program is run as follows:

```
java -cp talairach.jar org.talairach.PointToTD 1, 15, 10, 12
```

In this example, the Structural Probability Map results for (15,10,12) would be returned. The cube range search uses 5mm (or +/-2mm) as its default. To set a different cube size, use "3:<cube size>". Sizes of 3, 5, 7, 9 and 11 are accepted. For example:

```
java -cp talairach.jar org.talairach.PointToTD 3:7, 15, 10, 8
```

That search would return the nearby 343 labels (for a 7x7x7 region).

## ExcelToTD

ExcelToTD operates on a tab-delimited text file instead of a single point. It has the same search options as PointToTD. Unlike the graphical Talairach Client, ExcelToTD does not have the option of specifying a row and column where coordinate data starts. It assumes the first three items on each line are "x y z". For example:

```
java -cp talairach.jar org.talairach.ExcelToTD 2, data.txt
```

Each coordinate listed in data.txt will be assigned a Talairach label (option 2). The coordinates and labels will be saved in a tab-delimited file named data.td.

## Talairach Client

The Talairach Client can also be launched via the command line:

```
java -cp talairach.jar org.talairach.TalairachClient
```

## Suggested Uses

### Brodmann Area (BA) Labels for Cortical Activations

When seeking Brodmann Area labels for cortical activation sites it is possible to extend the search diameter to find the nearest gray matter labels. For an experiment designed to activate the M1 mouth motor region only 38% of the sites were found to fall within Brodmann Areas, but as the search diameter was increased to 3 mm, 5 mm, and 7 mm, BA labels were obtained for 62%, 92%, and 100% of the sites. This example shows the utility of obtaining BA labels for cortical activation sites. While the appropriate search range may vary from site to site within the brain, this is easy to test. As the search range increases the number of labels found increases. This presents two problems to the user. First, the label retrieval process becomes much slower. For a 5x5x5 mm search range 125 voxels are searched and many labels are found. The Talairach Client reduces the number of labels by only responding with unique labels. The unique labels are organized by incidence within the search range, with the highest label incidence being the first label returned. Along with each of the unique labels is the number of voxel within the search range with that label. Second, with large search ranges, the user has to deal with much larger files to sort through. The recommended strategy is to search with a small search region initially, remove coordinates with BA labels, and proceed using smaller input files for larger search ranges.

### Anatomical Organization of Coordinate Data

For applications that provide Talairach coordinates (FSL, SPM, AFNI, MEDx, TurboFire, etc.), it is helpful to use anatomical labels to organize findings anatomically. For applications that provide only MNI coordinates, users should utilize the [icbm2tal transform](#) to convert

coordinates between MNI space and Talairach space. The data saved by the Talairach Client can be used to create a labeled file. A common use of this labeled file is to open it in Excel and rearrange the original slice-ordered data into 3-D anatomical groupings by sorting by lobe and gyrus. This provides users with a good sense of which activation sites to group together anatomically.

Copyright © 2003-2015 Research Imaging Institute. All rights reserved.

Converting Between MNI and Talairach Spaces

In the past, BrainMap has utilized the Brett transform to convert MNI coordinates to Talairach space ([mni2tal](#)). However, we recently switched to a new MNI space to Talairach space transform called **icbm2tal**.

The findings of [Lancaster et al., 2007](#), [📄](#) show that MNI/Talairach coordinate bias associated with reference frame (position and orientation) and scale (brain size) can be substantially reduced using the best-fit icbm2tal transform. This transform has been validated and shown to provide improved fit over the Brett mni2tal transform.

Separate best-fit transforms were determined for SPM2 (**icbm\_spm2tal**) and FSL (**icbm\_fsl2tal**). A best-fit icbm2tal transform was also formulated pooling SPM2 and FSL data (**icbm\_other2tal**).

Implementation of icbm2tal in BrainMap

In BrainMap, we record both the brain template and the software used for spatial normalization. When this information is inserted into the database, we apply a coordinate transformation based on the reported template/software combination.

Currently, a transform is applied to all MNI coordinates that are added to the database. If the published coordinates were converted from MNI space to Talairach space using the Brett transform, then BrainMap will "unBrett" these coordinates using tal2mni, and apply a version of the icbm2tal prior to insertion into the database.

Reported Template

Brett Transform – FSL  
Brett Transform – SPM5  
Brett Transform – SPM2  
Brett Transform – SPM96  
Brett Transform – SPM97  
Brett Transform – SPM99  
MNI – AFNI  
MNI – FSL  
MNI – In-House  
MNI – SPM5  
MNI – SPM2  
MNI – SPM96  
MNI – SPM97  
MNI – SPM98  
MNI – SPM99  
MNI – Unknown SPM

Transformation Applied in BrainMap

Brett tal2mni, then icbm\_fsl2tal  
Brett tal2mni, then icbm\_spm2tal  
Brett tal2mni, then icbm\_spm2tal  
Brett tal2mni, then icbm\_spm2tal  
Brett tal2mni, then icbm\_spm2tal  
Brett tal2mni, then icbm\_spm2tal  
icbm\_other2tal  
icbm\_fsl2tal  
icbm\_other2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal  
icbm\_spm2tal

Have a question? Check our forums!

- [brainmap.org/forum](#)

Find a previous answer or ask the experts yourself

**Quick Author Search**

Want to check if a paper is already in the BrainMap database? Just type in the author's last name below:

Search

**Activation Coordinate  
Experiment-wise Search (ACES)**

Upload a tab-delimited file of locations to find which BrainMap experiments are most similar:

Choose File

No file chosen

Reference space: ☒ Talairach ☐ MNI

Find similar experiments: 

Search


**Functional Database Status**

Papers: 3026  
Experiments: 14887  
Paradigm Classes: 108  
Subjects: 64277  
Locations: 116639

**VBM Database Status**

Papers: 980  
Experiments: 3093  
Subjects: 73938  
Locations: 21481

**Current Software Versions**

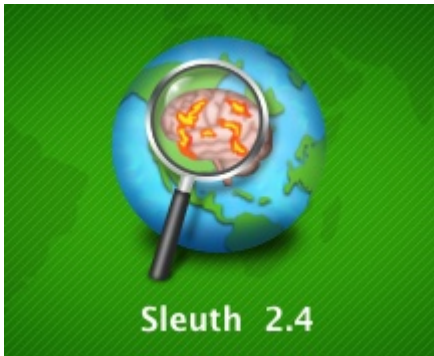

**PLEASE NOTE:** The transformations for converting between MNI coordinates obtained using SPM2 and FSL have been validated. All of the other ways in which we transform MNI coordinates to Talairach space (e.g., SPM5, SPM96, SPM99, etc.) have not been validated and therefore should be interpreted with caution. We have tried to make the best of what can only be viewed as a nasty situation.

In contrast to MNI coordinates, BrainMap does not perform any transformation of coordinates in Talairach space (listed below). We cannot be sure if this is the correct approach for all software packages. Specifically, we are sure that some transformation should be applied to coordinates derived from the Human Brain Atlas ([Roland and Zilles, 1994](#)). However, we are unable to do so until a transformation has been developed and validated.

#### Reported Template

Talairach 1988

Talairach 1988 – AFNI

Talairach 1988 – AIR

Talairach 1988 – Brain Voyager

Talairach 1988 – BRAINS

Talairach 1988 – BrainVOX

Talairach 1988 – Human Brain Atlas

Talairach 1988 – LIPSIA

Talairach 1988 – MedX

Talairach 1988 – SPM 4.0

Talairach 1988 – SPM94


Talairach 1988 – SPM95

Talairach 1993

**ONE LAST THING TO CONSIDER:** We have attempted to accurately record the template and software used for spatial normalization. However, this information can be difficult to ascertain. Specifically, when authors use SPM for data analysis and report their coordinates as "Talairach coordinates", it is often very hard to determine if they used the Brett transform or are simply unaware of the differences between MNI and Talairach spaces. We have done the best that we can. If you see an error in the recorded template for a given paper, please email us and we will correct our mistake.

#### icbm2tal Downloads

To convert your coordinates between the MNI and Talairach spaces, you can use our java application, [GingerALE](#) to convert your coordinates. In GingerALE, click on Tools → Convert Foci. This menu item will open a dialog window that will guide you through the conversion of your coordinates.

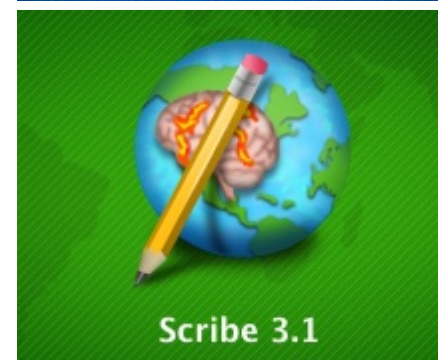





Alternatively, you can download the appropriate Matlab file:


- [icbm\\_spm2tal.m](#)
- [icbm\\_fsl2tal.m](#)
- [icbm\\_other2tal.m](#)
- [tal2icbm\\_spm.m](#)
- [tal2icbm\\_fsl.m](#)
- [tal2icbm\\_other.m](#)

Copyright © 2003-2016 Research Imaging Institute. All rights reserved.

[home](#) [client](#) [applet](#) [daemon](#) [about](#)

## Talairach Daemon

The Talairach Daemon (TD) is a high-speed database server for querying and retrieving data about human brain structure over the internet. It uses TCP/IP sockets for communications and minimizes the amount of data transferred during transactions. By keeping most transactions to a low number of bytes (less than 50 generally), even slow throughput network transfers (1 Kbyte/sec) should have reasonable response times.

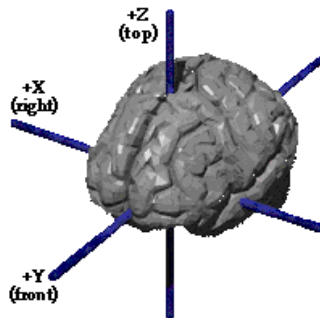

The TD server data is searched using x,y,z coordinates resolved to 1x1x1 mm volume elements within a standardized Talairach space. An array, indexed by x,y,z coordinates, that spans the dimensions of the 1988 Talairach Atlas brain (170mm, 210mm, 200mm), provides high-speed access to data. Coordinates tracked by the TD server are spatially consistent with the Talairach Atlas. Each array location stores a pointer to a relation record that holds data describing what is present at the corresponding coordinate. Presently, the data in relation records are either Structure Probability Maps (SP Maps) or Talairach Atlas Labels, though others can be easily added. The relation records are implemented as linked lists to names and values for brain structures.

The Daemon can be accessed using the [Talairach Applet](#) or by using the [Talairach Client](#).

To run a local copy of the Talairach Daemon, download the [.jar](#) file and run the following command:

```
java -cp talairach.jar org.talairach.AtlasServer 1600
```

To access a locally-run Daemon, use the [command line](#) options with an additional parameter: host=www.yourserver.com:port.

support from the [EJLB Foundation](#) and the [NIH Human Brain Project](#).

This work forms part of a continuing project of the [ICBM](#) to develop a probabilistic atlas of human neuroanatomy.

---

### Recommended System Requirements

- [Java](#) 1.4 or better
- 256 MB RAM

---

### Java Runtime Environment

Update the [Java Runtime Environment](#) (JRE) at Sun's official site.

---

### Feedback

Please email [Mick Fox](#) if you have any questions or find any bugs while using this software.